Laboratorní a zdravotnická technika OPTING servis Bohumínská 788/61 710 00 Ostrava 10

telefon: 596 241 785, 596 241 852-i fax, 596 241 960, 603 464 790 optingservis@optingservis.cz www.optingservis.cz

General Medicine and ENT Instruments

01

THE BEST IN QUALITY AND PERFORMANCE

High tech materials, the latest production processes and innovative design are three factors which contribute to the excellent quality of HEINE Otoscopes. They produce excellent image quality with uniquely bright and reflex-free illumination of the tympanum and auditory canal.

HEINE Otoscopes are available for various applications with distal fiber optic (F.O.) or direct bulb illumination with bright HEINE XHL® Xenon halogen bulbs.





HEINE Otoscopes



K 180

mini 3000 F.O.

BETA 100

K 100

Operating otoscope

mini 3000

Metal Polycarbonate

Glass Acrylic 3x magnification

> F.O. Direct

Professional Compact

Battery 2.5 V Rechargeable battery 2.5 V Rechargeable battery 3.5 V

Reusable tips Reusable Specula AllSpec Disposable tips UniSpec Disposable Specula Tip-Adaptor

> Insufflation port Page

		Hou	sing construc	tion		
		٧	iewing windo	w		
					2,5 x	
			Illumination			
		ı	nstrument siz	:e		
		Availa	ıble power so	urces		
			Tips			
			*	□*	*	
			sufflation Bul			
		□**				
009	011	012	015	017	019	020

^{*} with Tip-Adapter
** Optional insufflation port connector available



High-quality Fiber Optic (F.O.) illumination with high-intensity XHL Xenon Halogen bulb.







HEINE BETA 200® Fiber Optic (F.O.) Otoscope

The ultimate in performance and quality

Solid, all-metal construction. A quality instrument for the physician's office with Fiber Optic (F.O.) illumination.

- : All-metal housing matt-chrome-plated. Durable, long-lasting and easy to clean.
- : Instrument head matt-black inside. Eliminates reflexes.
- : Distal Fiber Optic (F.O.) illumination with XHL Xenon Halogen Technology. Homogenous, reflex-free illumination of ear canal and tympanum.
- : Large, scratch-resistant glass viewing window with 3x magnification. Unobstructed view, distortion-free image.
- : Airtight window seal. Reliable test of tympanic mobility.
- : Swivelling viewing window built into instrument. Useful for instrumentation, cannot be mislaid.

BETA200 F.O. Otoscope	2.5V	3.5 V
with 4 reusable tips, without handle	B-001.11.501	B-002.11.501
same, without tips, without handle	B-001.11.500	B-002.11.500
XHL Xenon Halogen spare bulb	X-001.88.077	X-002.88.078

For information on diagnostic combined sets, please see pages 157 – 159.



- : Fiber Optic (F.O.) illumination
- : All-metal housing
- : 3x magnification
- : Silica glass viewing window

HEINE BETA 200® Fiber Optic (F.O.) ENT Diagnostic Sets



spare bulb hard case



BETA 200 F.O. Otoscope Sets	2.5 V	3.5 V (NiMH)	3.5 V II ION
BETA battery handle	B-141.10.118		
BETA R rechargeable handle for			
mains socket		B-141.20.376	B-141.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		B-141.20.384	B-141.27.384
BETA NT rechargeable handle		-	
and NT 300 charger		B-141.20.420	
BETA L rechargeable handle			
and NT 300 charger			B-141.29.420

(III): Rechargeable handle with charge status display in the bottom insert.



HEINE BETA 200® Fiber Optic (F.O.) Diagnostic Sets



1 AllSpec set of reusable tips (B-000.11.111) 10 disposable tips and one spare bulb each hard case

BETA 200 Diagnostic Sets	2.5 V	3.5 V (NiMH)	3.5V LI ION
BETA battery handle	A-132.10.118		
BETA R rechargeable handle for			
mains socket		A-132.20.376	A-132.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		A-132.20.384	A-132.27.384
BETA NT rechargeable handle			
and NT 300 charger		A-132.20.420	
BETA L rechargeable handle			
and NT 300 charger			A-132.29.420

As above, but with additional multi-purpose illuminator



BETA 200 Diagnostic Sets	2.5V	3.5V (NiMH)	3.5V LI ION
BETA battery handle	A-134.10.118		
BETA R rechargeable handle for			
mains socket		A-134.20.376	A-134.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		A-134.20.384	A-134.27.384
BETA NT rechargeable handle			
and NT 300 charger		A-134.20.420	
BETA L rechargeable handle			
and NT 300 charger			A-134.29.420

①: Rechargeable handle with charge status display in the bottom insert.



HEINE K 180[®] Fiber Optic (F.O.) Otoscope

The standard, all around Fiber Optic Otoscope



Otoscope with Fiber Optic (F.O.) illumination constructed of light, high-impact polycarbonate and $3\,x$ magnification scratch-resistant acrylic viewing window. Otherwise technically identical to the BETA200 F.O. Otoscope.

K180 F.O. Otoscope	2.5 V	3.5 V
with 4 reusable tips, without handle		B-002.11.551
same, without tips, without handle	B-001.11.550	B-002.11.550
XHL Xenon Halogen spare bulb	X-001.88.077	X-002.88.078

: Fiber Optic (F.O.) illumination

: Polycarbonate construction

HEINE K 180® Fiber Optic (F.O.) Diagnostic Sets



Set complete with: **K180 Fiber Optic (F.O.) Otoscope** 1 set of reusable tips (B-000.11.111) soft pouch

K180 F.O. Otoscope Sets	2.5 V	3.5 V (NiMH)	3.5V LI ION
BETA battery handle	B-181.10.118		
BETA R rechargeable handle for			
mains socket		B-181.20.376	B-181.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		B-181.20.384	B-181.27.384
BETA NT rechargeable handle			
and NT 300 charger		B-181.20.420	
BETA L rechargeable handle			
and NT 300 charger			B-181.29.420



Set complete with K180 Ophthalmoscope, K180 Fiber Optic (F.O.) Otoscope 1 set of reusable tips (B-000.11.111) hard case

K180 Combined Diagnostic Sets	2.5 V	3.5V (NiMH)	3.5V LI ION
BETA battery handle	A-279.10.118		
BETA R rechargeable handle for			
mains socket		A-279.20.376	A-279.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		A-279.20.384	A-279.27.384
BETA NT rechargeable handle			
and NT 300 charger		A-279.20.420	
BETA L rechargeable handle			
and NT 300 charger			A-279.29.420



HEINE mini 3000® Fiber Optic (F.O.) Otoscope

Compact pocket otoscope with Fiber Optic (F.O.) illumination

Modern pocket otoscope in a unique, compact design. Handy and powerful. Can only be used with the mini 3000 handle system. Available in black or blue.

- : Instrument head matt-black inside. Eliminates reflexes.
- Distal Fiber-Optic (F.O.) illumination with enhanced XHL Xenon Halogen Technology producing 100% more light compared with conventional bulbs. Homogenous, reflex-free illumination of ear canal and tympanum.
- : 2-piece handle and head. Easy to maintain, flexible. Compatible with other mini 3000 instrument heads.
- : Insufflation port. For pneumatic testing of tympanic mobility. Please order connector separately.
- : Viewing window with 3x magnification and optimized casing surface for razor-sharp images and minimal reflection. Swivelling viewing window built into instrument.
- : High-quality handle: Chrome-finish upper section/refined plastic. Shockproof, sturdy, non-slip.
- : Attachment clip with integrated on/off switch. Secure. Switches off automatically when returned to the pocket. 20.000 switch cycles guaranteed.
- : Replaceable batteries. Size AA, or optional rechargeable battery with mini NT charger*.

mini 3000 F.O. Otoscope	2.5 V
mini 3000 F.O. Otoscope with mini 3000 battery handle and 1 set	
(4 pcs.) reusable tips and 5 of both 2.5 and 4 mm disposable tips, with	
batteries	D-001.70.120
mini 3000 F.O. Otoscope with mini 3000 battery handle and 5 of both	
2.5 and 4 mm disposable tips, with batteries	D-001.70.110
mini 3000 F.O. Otoscope with bulb, without handle and without	
accessories	D-001.70.106
XHL Xenon Halogen spare bulb	X-001.88.105
Insufflation port connector for mini 3000 F.O.	D-000.80.101



: Modern pocket otoscope

: Fiber Optic (F.O.) illumination

: Switches off automatically

: Optional rechargeable handle



*To upgrade to a rechargeable instrument, with mini NT charger, please see page 142.



HEINE mini 3000® pocket instruments – a complete range of quality diagnostic instruments for your pocket. Available in black or blue. Please specify the colour you prefer when ordering. If not otherwise specified, we will supply black instruments.



HEINE mini 3000® Diagnostic Sets



Set D-851 complete with: mini 3000 Fiber Optic (F.O.) Otoscope 1 set = 4 reusable tips (B-000.11.111) 10 AllSpec disposable tips hard case

mini 3000 Diagnostic Sets	2.5 V
mini 3000 battery handle with batteries	D-851.10.021
as above, but with mini 3000 Otoscope	D-851.20.021



Set D-873 complete with: mini3000 Ophthalmoscope, mini3000 F.O. Otoscope, 1 set = 4 reusable tips (B-000.11.111)
5 each of 2.5 and 4 mm dia. AllSpec disposable tips hard case

mini 3000 Diagnostic Sets	2.5 V
2 mini 3000 battery handles with batteries	D-873.11.021
same but with mini 3000 Otoscope	D-873.21.021



*To upgrade to a rechargeable instrument, with mini NT charger, please see page 142.

HEINE AllSpec® Disposable tips







- \div For single use. Hygienic, no cross-infection.
- \div Original HEINE design. Firm fit to the otoscope head.
- : Assured quality. No sharp edges.

AllSpec Disposable tips (grey)	
for BETA 200, K180, mini 3000 F.O., mini 3000 Otoscopes	
Pack of 1000 pcs. 2.5 mm dia. (infants) [03]	B-000.11.128
Pack of 1000 pcs. 4mm dia. (adults) [03]	B-000.11.127
Case of 10 dispenser packs of 250 tips, 2.5 mm dia (infants) [01]	B-000.11.151
Case of 10 dispenser packs of 250 tips, 4 mm dia (adults) [01]	B-000.11.150
AllSpec Tip-dispenser [02]	B-000.11.148

AllSpec Tip-dispenser [02]

The right size, visible and instantly available. A clean new tip for each and every patient. Contents: $50 \text{ pcs. } \emptyset 2.5 \text{ mm}$ and $50 \text{ pcs. } \emptyset 4 \text{ mm}$.

Ø 10 mm



Reusable tips



- : Made of high-density SANALON S plastic with smooth external finish. Easy to clean. All common germicides can be used. Autoclavable at up to 134°C.
- : Ergonomic shape. Painless atraumatic insertion.
- : Positive locking design. Firm fit to the otoscope head.



Soft Disposable tips for reusable tips

: Fits over the end of the AllSpec tip to protect the auditory canal.



Soft disposable tips (black)	
Soft disposable tip 3 mm dia. 40 pcs./1 pack	B-000.11.141
Soft disposable tip 5 mm dia. 40 pcs./1 pack	B-000.11.142

Insufflation Bulb




: For pneumatic testing.

Insufflation bulb for BETA 200 , K180, BETA 100 and K100 Otoscopes	B-000.11.240
Insufflation bulb for mini 3000 F.O. Otoscope with connector	D-000.80.102
Insufflation port connector for mini 3000 F.O. Otoscope	D-000.80.101



Bright XHL Xenon Halogen direct illumination.




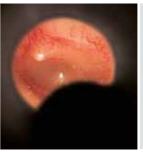



HEINE BETA 100® Diagnostic Otoscope

Diagnostic Otoscope with an all-metal housing

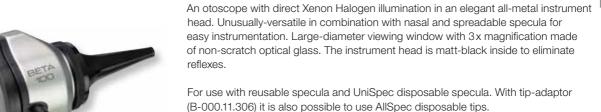

- : All-metal housing. Exceptionally tough and durable.
- : Instrument head matt-black inside. Eliminates reflexes.
- : XHL Xenon Halogen Technology. Extremely bright, white light.
- : Large, scratch-resistant glass viewing window with 3x magnification. Unobstructed view, distortion-free image.
- : Swivelling viewing window built into instrument. Useful for instrumentation, cannot be mislaid.
- : Insufflation port. Reliable test of tympanic mobility.
- : Extensive choice of accessories. Very versatile e. g. for nasal examination.

BETA100 Diagnostic Otoscope	2.5V	3.5 V
with 6 reusable specula, without handle without specula, without handle		B-002.11.527 B-002.11.525
XHL Xenon Halogen spare bulb	X-001.88.037	X-002.88.049
Replacement sleeve for bulb. Reduces reflexes. [01]		B-000.11.212



- : All-metal housing
- : Silica glass viewing window
- : Easy instrumentation





HEINE BETA 100® ENT Diagnostic Sets



Set complete with: **BETA100 Diagnostic Otoscope** 1 set (6 pcs.) of reusable specula (B-000.11.221) 10 UniSpec disposable tips spare bulb hard case

BETA 100 Otoscope Sets	2.5V	3.5 V (NiMH)	3.5V LI ION
BETA battery handle	B-236.10.118		
BETA R rechargeable handle for			
mains socket		B-236.20.376	B-236.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		B-236.20.384	B-236.27.384
BETA NT rechargeable handle			
and NT 300 charger		B-236.20.420	
BETA L rechargeable handle			
and NT 300 charger			B-236.29.420

①: Rechargeable handle with charge status display in the bottom insert.

HEINE BETA Diagnostic Set

Set complete: **BETA 200 Ophthalmoscope, BETA 100 Diagnostic Otoscope,** 1 set (6 pcs.) of reusable specula (B-000.11.221) 10 UniSpec disposable tips one spare bulb each hard case



BETA Diagnostic Sets	2.5V	3.5V (NiMH)	3.5V LI ION
BETA battery handle	A-234.10.118		
BETA R rechargeable handle for			
mains socket		A-234.20.376	A-234.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		A-234.20.384	A-234.27.384
BETA NT rechargeable handle			
and NT 300 charger		A-234.20.420	
BETA L rechargeable handle			
and NT 300 charger			A-234.29.420

(I): Rechargeable handle with charge status display in the bottom insert.



HEINE K 100® Diagnostic Otoscope

Diagnostic Otoscope with a high-impact polycarbonate housing



A very cost-effective diagnostic otoscope with direct Xenon Halogen illumination. The head is moulded in high-impact polycarbonate, and features a high quality acrylic viewing window with 3x magnification. Matt black interior to eliminate reflexes.

For use with reusable specula and UniSpec disposable specula. With tip-adaptor (B-000.11.306) it is also possible to use AllSpec disposable tips.

- : Polycarbonate construction. Shock-resistant.
- : Instrument head matt-black inside. Eliminates reflexes.
- : XHL Xenon Halogen Technology. Extremely bright, white light.
- : Large acrylic viewing window with **3x magnification**. Sharply-focused distorsion-free image.
- : Swivelling viewing window built into instrument. Useful for instrumentation, cannot be mislaid.
- : Insufflation port. Reliable test of tympanic mobility.

: Polycarbonate construction
: Easy instrumentation
: Easy instrument

K100 Diagnostic Otoscope	2.5 V	3.5 V
with 3 reusable specula 3.5/4.5/5.5mm dia., without handle same, without specula, without handle		B-002.11.576 B-002.11.575
XHL Xenon Halogen spare bulb	X-001.88.037	X-002.88.049
Replacement cover for bulb [01] Reduces reflexes.		B-000.11.212

HEINE K 100® ENT Diagnostic Sets

Set complete with: **K100 Diagnostic Otoscope** 1 set (6 pcs.) of reusable specula (B-000.11.221) spare bulb soft pouch



Sets K100 Otoscope	2.5V	3.5 V (NiMH)	3.5V LI ION
BETA battery handle	B-237.10.118		
BETA R rechargeable handle for			
mains socket		B-237.20.376	B-237.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		B-237.20.384	B-237.27.384
BETA NT rechargeable handle			
and NT 300 charger		B-237.20.420	
BETA L rechargeable handle			
and NT 300 charger			B-237.29.420

(II): Rechargeable handle with charge status display in the bottom insert.



HEINE K 100® ENT Diagnostic Sets

Set complete with: **K100 Diagnostic Otoscope** 3 reusable specula (B-000.11.222) soft pouch




Sets K100 Otoscope	2.5V	3.5V (NiMH)	3.5V IIIN
BETA battery handle	B-245.10.118		
BETA R rechargeable handle for			
mains socket		B-245.20.376	B-245.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		B-245.20.384	B-245.27.384
BETA NT rechargeable handle			
and NT 300 charger		B-245.20.420	
BETA L rechargeable handle			
and NT 300 charger			B-245.29.420

: Rechargeable handle with charge status display in the bottom insert.



HEINE Operating Otoscope

The greatest flexibility and most options for instrumentation

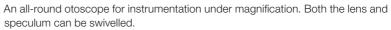

Ideal for use with nasal and spreadable specula.

For use with reusable specula and UniSpec disposable specula. With tip-adaptor (B-000.11.306) it is also possible to use AllSpec disposable tips.

- : Open construction. Ideal for instrumentation.
- : XHL Xenon Halogen Technology. Extremely bright, white light.
- : Swivelling lens and speculum. Clears the way for access and examination under magnification and illumination.
- : 2.5x lens in protective frame. Prevents scratching.
- : The relative positions of light bulb and speculum can be fixed. Helps with instrumentation, easily returned to normal position.

HEINE Operating Otoscope	2.5 V	3.5 V
with 1 set = 6 reusable specula, without handle	B-001.11.494	B-002.11.494
same, without specula, without handle	B-001.11.492	B-002.11.492
XHL Xenon Halogen spare bulb	X-001.88.037	X-002.88.049

HEINE Operating Otoscope ENT Diagnostic Sets



Set complete with: **Operating Otoscope**1 set (5 pcs.) of reusable specula (B-000.11.221)
10 UniSpec disposable tips
spare bulb
hard case

Operating Otoscope Sets	2.5 V	3.5 V (NiMH)	3.5V LI ION
BETA battery handle	B-188.10.118		
BETA R rechargeable handle for			
mains socket		B-188.20.376	B-188.27.376
BETA TR rechargeable handle			
and transformer with plug-in			
transformer		B-188.20.384	B-188.27.384
BETA NT rechargeable handle			
and NT 300 charger		B-188.20.420	
BETA L rechargeable handle			
and NT 300 charger			B-188.29.420

①: Rechargeable handle with charge status display in the bottom insert.

: Modern pocket otoscope

Switches off automaticallyOptional rechargeable handle



HEINE mini 3000® Otoscope

Compact pocket otoscope with direct illumination

Modern pocket otoscope in a unique, compact design. Can only be used with the mini 3000 handle system.

Available in black or blue.



- : Enhanced XHL Xenon Halogen Technology with 100% more light compared with conventional bulbs. Very bright, concentrated white light for perfect illumination.
- **2-piece handle and head.** Easy to maintain, flexible. Compatible with other mini 3000 instrument heads.
- : Viewing window with 3x magnification and optimized casing surface for razor-sharp images and minimal reflection.
- : Swivelling viewing window built into instrument. Useful for instrumentation, cannot be mislaid.
- : High-quality handle: Chrome-finish upper section/refined plastic. Shockproof, sturdy, non-slip.
- : Attachment clip with integrated on/off switch. Secure. Switches off automatically when returned to the pocket. 20.000 switch cycles guaranteed.
- : Replaceable batteries. Size AA, or optional rechargeable battery with mini NT charger*.

mini 3000 Otoscope	2.5 V
mini 3000 Otoscope with mini 3000 battery handle and 1 set (4 pcs.)	
reusable tips and 5 of both AllSpec disposable tips 2.5 and 4 mm dia.,	
with batteries	D-001.70.220
mini 3000 Otoscope with mini 3000 battery handle and 5 of both	
AllSpec disposable tips 2.5 and 4 mm dia., with batteries	D-001.70.210
mini 3000 Otoscope with bulb, without handle and without	
accessories	D-001.70.206
XHL Xenon Halogen spare bulb	X-001.88.110

For information on combined sets, please see pages 157 – 159.




*To upgrade to a rechargeable instrument, with mini NT charger, please see page 142.



HEINE mini 3000® pocket instruments – a complete range of quality diagnostic instruments for your pocket. Available in black or blue. Please specify the colour you prefer when ordering. If not otherwise specified, we will supply black instruments.



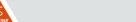

HEINE UniSpec® Disposable Ear Specula



- : Disposable. Hygienic and safe.
- : Original HEINE design and manufacture. Fit positively to the otoscope.
- : Assured quality. No sharp edges.
- : A clean, new speculum for each and every patient. The right size visible and instantly available.

UniSpec Disposable Ear Specula (grey)	
for BETA 100, K100, Operating Otoscope	
Pack of 1000, 2.5 mm dia. (infants) 1000 pcs./1 pack	B-000.11.242
Pack of 1000, 4 mm dia. (adults) 1000 pcs./1 pack	B-000.11.241
UniSpec specula dispenser	B-000.11.146

Reusable Specula





4.5 mm dia.







- : Moulded in high-density SANALON S, smooth surface. Easy to clean, can be disinfected with conventional solutions. Autoclavable at 134°C.
- : Ergonomic Design. Painless and atraumatic insertion.
- : Metal security lock with bayonet fitting. Long-lasting, positive fit to the otoscope.

Reusable Specula (black)	
for BETA100, K100 and Operating Otoscope	
1 set = 5 ear specula and 1 nose speculum	B-000.11.221
3 specula 3.5/4.5/5.5mm dia.	B-000.11.222
Ear speculum 2.2 mm dia.	B-000.11.215
Ear speculum 2.8 mm dia.	B-000.11.216
Ear speculum 3.5 mm dia.	B-000.11.217
Ear speculum 4.5 mm dia.	B-000.11.218
Ear speculum 5.5 mm dia.	B-000.11.219
Nose speculum 10 mm dia.	B-000.11.220

Spreadable Nasal Speculum



- : All-metal construction. Anatomically-correct rounded shape. Non-injuring.
- : For use with Operating Otoscope, BETA 100 Diagnostic Otoscope and Fiber Optic Nasal Illuminator.

Spreadable Nasal Speculum	B-000.11.231
---------------------------	--------------

Not recommended for use with K100 Diagnostic Otoscopes.



Universal Speculum



- : All-metal construction. Opens symmetrically from 2.5 to 15 mm.
- : Rotate ring to open branches. Stays put in any position.
- For use with Operating Otoscope, BETA100 Diagnostic Otoscope and Fiber Optic Nasal Illuminator.

Universal Speculum

B-000.11.239

Not recommended for use with K100 Diagnostic Otoscopes.

Tip-Adaptor



: All-metal construction. Adaptor for AllSpec disposable tips (p. 13) for use with diagnostic (BETA 100, K 100) and operating otoscope.

Tip-Adaptor B-000.11.306

Insufflation Bulb



For Diagnostic Otoscopes.

: For pneumatic testing.

Insufflation bulb for BETA 100 and K 100 Otoscopes

B-000.11.240



A company makes history and gives a name to outstanding product quality: HEINE.

Founded in 1946, family owned, rich in tradition and today international market leader in a complete range of diagnostic instruments. HEINE products meet all international standards (ISO/CE) and represent the leading edge in precision and ergonomic design. Ongoing dedication to research and development has created a solid and diverse technology base for continued quality and product leadership. Over 500 employees worldwide contribute to this success.

www.heine.com

HEINE Optotechnik GmbH & Co. KG 82211 Herrsching Kientalstr. 7 Tel. +49 (0) 81 52-38 0 Fax +49 (0) 81 52-3 82 02 E-Mail: info@heine.com

HEINE USA LTD. 10 Innovation Way Dover, NH 03820 Tel. (603) 7 42-71 03 Fax (603) 7 42-72 17 Toll Free (800) 367-4872 E-Mail: service@heine-na.com

HEINE AUSTRALIA PTY. LTD. Unit 9, 98 Old Pittwater Road P.O. Box 717 Brookvale 2100 New South Wales Tel. +61 (0) 2-99 38 95 00 Fax +61 (0) 2-99 39 23 05 E-Mail: info@heine.com.au

SWITZERLAND

HEINE (Schweiz) AG Tobeläckerstr. 9 CH-8212 Neuhausen am Rheinfall Tel. +41 (0) 52-6 72 22 66 Fax +41 (0) 52-6 72 63 77 E-Mail: heine-schweiz@schaffhausen.ch

We reserve the right to change

specification without notice.

Your distributor 10/10. A-000.00.100 e **Laboratorní a zdravotnická technika OPTING servis** Bohumínská 788/61 710 00 Ostrava 10

> telefon: 596 241 785, 596 241 852-i fax, 596 241 960, 603 464 790 optingservis@optingservis.cz www.optingservis.cz



